

legflow[®]

RX/OTW

Paclitaxel Releasing
Peripheral Balloon
Dilatation Catheter

Specifically designed for
SFA & BTK complex lesions






Coating technology

ReFlow Study:

International multicenter study in long (>20 cm) and complex lesions (TASC C & D)¹.

Legflow®'s performances in the treatment of long (>20 cm) TASC C&D lesions

- Primary Patency: 75.8% at 12 months and 70.5% at 24 months;
- Freedom from TLR: 82.2% at 12 months and 76.1% at 24 months;
- Survival rate: 94.7% at 12 months and 89.8% at 24 months.²

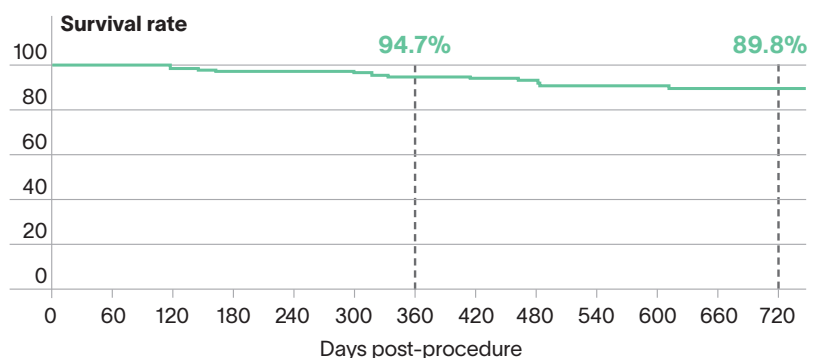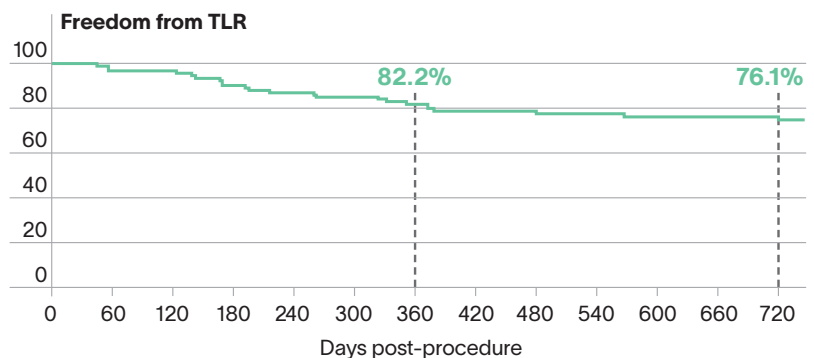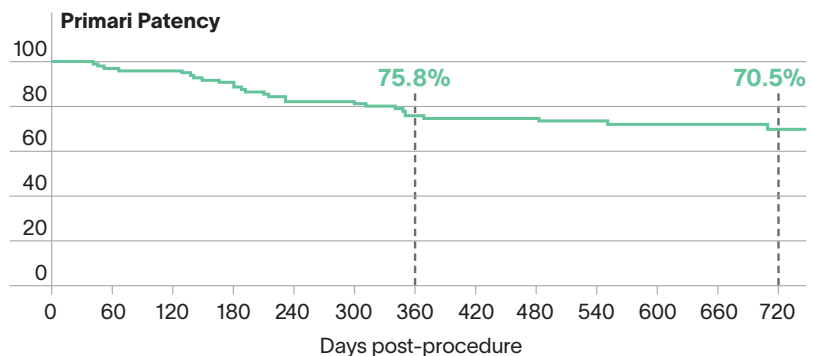


LegDEB Study:

International multicenter registry in real life population³

Legflow®'s advantages in complex patients with femoropopliteal lesions

Legflow® 1 year results in real life patients

Primary Patency: 81.1% and Freedom from TLR of 83.3% at 12 months;

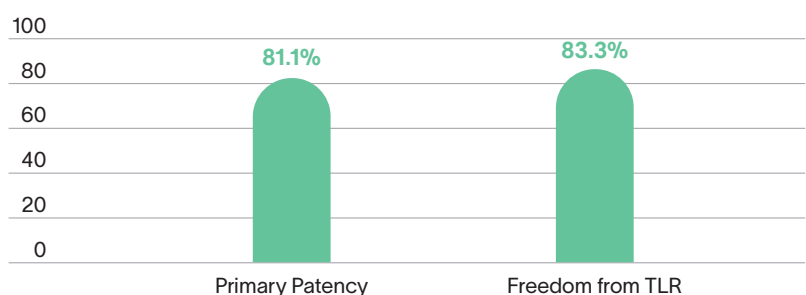

¹ Dr. Bosiers, Two-year outcome of the paclitaxel-eluting Legflow balloon catheter in the treatment of long and complex femoropopliteal lesions. The Journal of Cardiovascular surgery 2024 Dec 10 DOI: 10.23736/S0021-9509.24.12920-5


² Dr. Bosiers data presented @ Charing Cross 2021

³ Stabile E. et al., Journal of Endovascular Therapy 2019, Vol.26(1) 26-30.

Technical specifications

Technical data

A total of 163 options with:

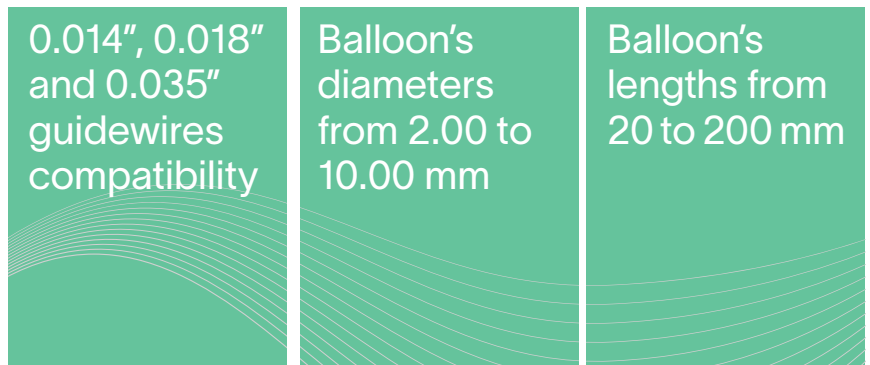

Type	OTW – 0.014"	OTW – 0.018"	OTW – 0.035"	RX – 0.014"
Shaft material	Polyamide	Polyamide	Polyamide	Stainless Steel
Balloon Material	Polyamide	Polyamide	Polyamide	Polyamide
Usable catheter length	150 cm	150 cm	80 or 135 cm	140 cm
Max. Recomm. Guidewire	0.014"	0.018"	0.035"	0.014"

Nominal Pressure	OTW – 0.014"	6 bar	
	OTW – 0.018"	6 bar	with lengths from 100 to 150 mm;
		8 bar	with lengths ≤ 80 mm
	OTW – 0.035"	6 bar	
Introducer sheath compatibility	RX – 0.014"	6 bar	
	OTW – 0.014"	4F	
	OTW – 0.018"	4F	Ø 2.0 to 4.0 mm with all lengths;
		5F	Ø 5.0 to 7.0 mm with all lengths.
	OTW – 0.035"	6F	Ø 4.00 to 6.00 mm
		7F	Ø 7.00 to 10.00 mm
	RX – 0.014"	5F	

Drug coating technology

Drug	Paclitaxel
Drug dose	3.0 µg / mm²
Delivery matrix	SAFEPAx®
Coated area	Cylindrical section of the balloon, exceeding the proximal and distal markers

Legflow[®] OTW
0.014"

Legflow[®] OTW
0.018"

Legflow[®] OTW
0.035"

Legflow[®] RX

Ordering
informations

Balloon length (mm)	Balloon Ø (mm) 0.014" OTW			
	Catheter length 150 cm			
	2.0 mm	2.5 mm	3.0 mm	3.5 mm
40	L 2.0-40 OTW	L 2.5-40 OTW	L 3.0-40 OTW	L 3.5-40 OTW
80	L 2.0-80 OTW	L 2.5-80 OTW	L 3.0-80 OTW	L 3.5-80 OTW
120	L 2.0-120 OTW	L 2.5-120 OTW	L 3.0-120 OTW	L 3.5-120 OTW
150	L 2.0-150 OTW	L 2.5-150 OTW	L 3.0-150 OTW	L 3.5-150 OTW

Balloon length (mm)	Balloon Ø (mm) 0.018" OTW					
	Catheter length 150 cm					
	2.0 mm	3.0 mm	4.0 mm	5.0 mm	6.0 mm	7.0 mm
40	L 18 2.0-40	L 18 3.0-40	L 18 4.0-40	L 18 5.0-40	L 18 6.0-40	L 18 7.0-40
60	L 18 2.0-60	L 18 3.0-60	L 18 4.0-60	L 18 5.0-60	L 18 6.0-60	-
80	L 18 2.0-80	L 18 3.0-80	L 18 4.0-80	L 18 5.0-80	L 18 6.0-80	-
100	L 18 2.0-100	L 18 3.0-100	L 18 4.0-100	L 18 5.0-100	L 18 6.0-100	-
120	L 18 2.0-120	L 18 3.0-120	L 18 4.0-120	L 18 5.0-120	L 18 6.0-120	-
150	L 18 2.0-150	L 18 3.0-150	L 18 4.0-150	L 18 5.0-150	L 18 6.0-150	-

Balloon length (mm)	Balloon Ø (mm) 0.035" OTW						
	Catheter length 80 cm (LS) or 135 cm (L)						
	4.0 mm	5.0 mm	6.0 mm	7.0 mm	8.0 mm	9.0 mm	10.0 mm
20	L 4.0-20	L 5.0-20	L 6.0-20	L 7.0-20	L 8.0-20	L 9.0-20	L 10.0-20
	LS 4.0-20	LS 5.0-20	LS 6.0-20	LS 7.0-20	LS 8.0-20	LS 9.0-20	LS 10.0-20
40	L 4.0-40	L 5.0-40	L 6.0-40	L 7.0-40	L 8.0-40	L 9.0-40	L 10.0-40
	LS 4.0-40	LS 5.0-40	LS 6.0-40	LS 7.0-40	LS 8.0-40	LS 9.0-40	LS 10.0-40
60	L 4.0-60	L 5.0-60	L 6.0-60	L 7.0-60	L 8.0-60	L 9.0-60	L 10.0-60
	LS 4.0-60	LS 5.0-60	LS 6.0-60	LS 7.0-60	LS 8.0-60	LS 9.0-60	LS 10.0-60
80	L 4.0-80	L 5.0-80	L 6.0-80	L 7.0-80	L 8.0-80	-	-
	LS 4.0-80	LS 5.0-80	LS 6.0-80	LS 7.0-80	LS 8.0-80	-	-
100	L 4.0-100	L 5.0-100	L 6.0-100	L 7.0-100	L 8.0-100	-	-
	LS 4.0-100	LS 5.0-100	LS 6.0-100	LS 7.0-100	LS 8.0-100	-	-
120	L 4.0-120	L 5.0-120	L 6.0-120	L 7.0-120	L 8.0-120	-	-
	LS 4.0-120	LS 5.0-120	LS 6.0-120	LS 7.0-120	LS 8.0-120	-	-
150	L 4.0-150	L 5.0-150	L 6.0-150	L 7.0-150	L 8.0-150	-	-
	LS 4.0-150	LS 5.0-150	LS 6.0-150	LS 7.0-150	LS 8.0-150	-	-

Balloon length (mm)	Balloon Ø (mm) 0.014" RX				
	Catheter length 140 cm				
	2.0 mm	2.5 mm	3.0 mm	3.5 mm	4.0 mm
20	L 2.0-20 RX	L 2.5-20 RX	L 3.0-20 RX	L 3.5-20 RX	L 4.0-20 RX
40	L 2.0-40 RX	L 2.5-40 RX	L 3.0-40 RX	L 3.5-40 RX	L 4.0-40 RX
60	L 2.0-60 RX	L 2.5-60 RX	L 3.0-60 RX	L 3.5-60 RX	-
80	L 2.0-80 RX	L 2.5-80 RX	L 3.0-80 RX	L 3.5-80 RX	L 4.0-80 RX
100	L 2.0-100 RX	L 2.5-100 RX	L 3.0-100 RX	L 3.5-100 RX	-
120	L 2.0-120 RX	L 2.5-120 RX	L 3.0-120 RX	L 3.5-120 RX	L 4.0-120 RX
150	L 2.0-150 RX	L 2.5-150 RX	L 3.0-150 RX	L 3.5-150 RX	-
200	L 2.0-200 RX	L 2.5-200 RX	-	-	-

Innovative Technology.
Proven Performance.

Not All DCBs Are The Same



Amorphous vs
Crystalline

The 3rd generation, unique amorphous paclitaxel matrix system with the highest coating stability on the market.

- **No Distal Embolization**

Designed with an amorphous and no-crystalline coating technology to prevent distal particles embolization ensuring safer outcomes, especially in CLTI patients.

- **No Slow Flow (SCF)**

Effectively eliminates the risk of Slow Flow and No-Reflow phenomena, promoting uninterrupted perfusion.

- **Unmatched Drug Retention**

Features the lowest drug wash-off rate and particle release among current DCBs, guaranteeing precise and sustained drug delivery directly to the target lesion.

- **Clinically Proven Performance**

Exceptional safety and efficacy demonstrated in thousands of patients, supported by robust randomized controlled trials (RCTs) with follow-up extending up to 5 years.



Comparison between Safepax®'s innovative no-crystalline coating (top) and traditional crystalline coatings used by competitors*

* Cardionovum® data on file



Manufacturer, Headquarter and International Sales Office

CARDIONOVUM GmbH
Am Bonner Bogen 2
53227 Bonn – Germany
phone +49 (0) 228-90 90 59-0
fax +49 (0) 228-90 90 59-20

CARDIONOVUM Schweiz
GmbH c/o BDO AG
Feldmoosstrasse 12
8853 Laachen – Switzerland
phone +41 79 560 91 93

E-MAIL info@cardionovum.com
WEB www.cardionovum.com



L.ML.01.01.25